



## Heavy Duty Electric Coin Cell Crimping Equipment Coin Cell Assembly Machine

Our Product Introduction

for more products please visit us on [lithiumbatterymachine.com](http://lithiumbatterymachine.com)

### Basic Information

- Place of Origin: China
- Brand Name: MRBEST
- Certification: CE
- Model Number: MR-DF160
- Minimum Order Quantity: 1 Set
- Price: Negotiable
- Packaging Details: Standard Export Wooden Packing
- Delivery Time: 7-20 working days
- Payment Terms: T/T, Western Union, MoneyGram
- Supply Ability: 500 sets per month



### Product Specification

- Name: Electric Coin Cell Crimper
- Sealing Mold: CR20 Crimping Mold (CR2016, CR2025, CR2032)
- Sealing Pressure: Max. 1200Kg, Digital Display Can Be Set
- Sealing Stroke: 20mm
- Structure: Heavy Duty Steel Platform
- Dimension: L250\*W185\*H450 Mm
- Highlight: Electric Coin Cell Crimping Equipment, Heavy Duty Coin Cell Crimping Equipment, Electric Coin Cell Assembly Machine



### More Images



## Product Description

Model	Electric Crimping and Disassemble Machine TOB-DF-160
Voltage	110V/220V AC±10% 50Hz/60Hz
Power	100W
Pressure	Max 1200KG, digital display can be set
Work stroke	20mm
Die	Standard die for CR20 series crimping, other dies is optional
Advantages	Precision controlled sealing pressure; the pressure is shown on the screen and it can be set freely.
Warranty	One year with lifetime support
Dimension(L*W*H)	250mm*185mm*450mm
Net weight	30KG



### Coin Cell Cripmer Electric Crimping Machine for Coin Cell Assembly

This electric coin cell crimping machine mainly use for laboratory battery material research, it can be also used for small scale button battery production in factory. It is equipped with electric driver, compare with hydraulic type, it is labor-saving. You can also change the die for dry powder tablets, wet powder tablets, press molding, riveting and other operations.

#### Features:


- \* Touch screen operation, easy to set parameters, high degree of automation.
- \* Set and display the sealing pressure value is available.
- \* It can be put into glove box through  $\phi 360$  antechamber without barrier.
- \* The use of electric drive, coin cell cases crimping process without vibration and leakage.
- \* Precise mold ensures accurate and reliable sealing, strong steel structure design, stable and safe.
- \* Small size, easy operation, it can be operated in the glove box.
- \* Beautiful appearance, easy maintenance, It has the advantages of small size, convenient operation, accurate molding and so on.

#### Specifications

Model	Electric Crimping and Disassemble Machine TOB-DF-160
Voltage	110V/220V AC±10% 50Hz/60Hz
Power	100W
Pressure	Max 1200KG, digital display can be set
Work stroke	20mm
Die	Standard die for CR20 series crimping, other dies is optional
Advantages	Precision controlled sealing pressure; the pressure is shown on the screen and it can be set freely.
Warranty	One year with lifetime support
Dimension(L*W*H)	250mm*185mm*450mm
Net weight	30KG

### Picture Display



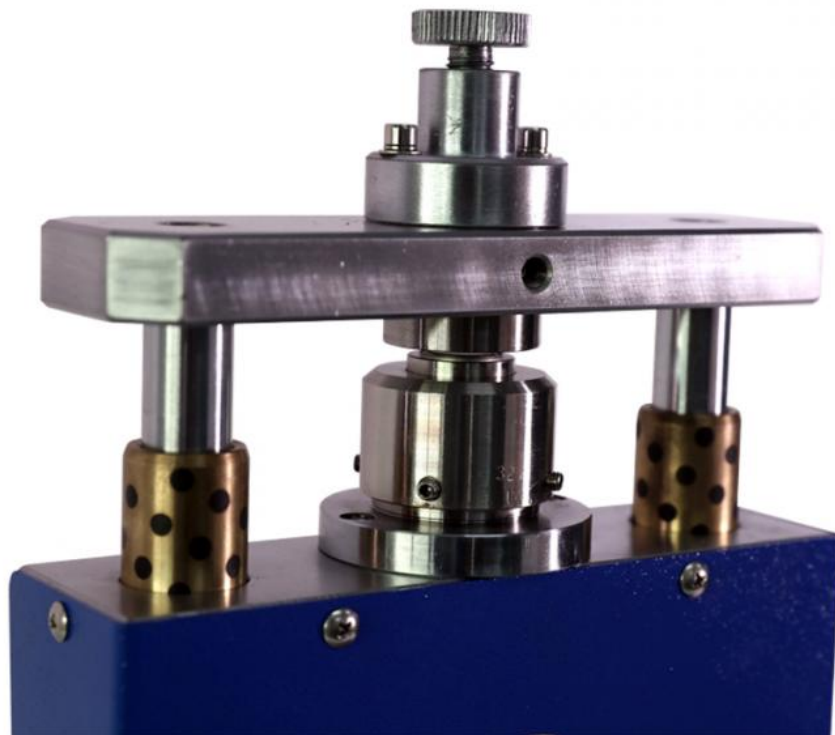













Planetary Vacuum Mixer

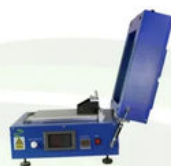

Automatic Film Coating Machine

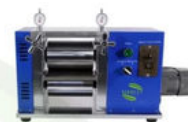

Rolling Press Machine



Coin Cell Punching Machine

### Coin Cell Lab Equipment



Vacuum Drying Oven



Glove Box



Coin Cell Crimper



Battery Testing Machine

## Factory Show

# COMPANY PROFILE

## 20 Years Lithium Battery Industry Experience



Shenzhen Meirui Zhida Technology Co., Ltd. is a Manufacturer specializing in a Complete Set of Lithium Ion Battery Research Equipment. MRBEST has been committed to providing customers with comprehensive Battery Research and Development Equipment to meet customers' new Lithium Battery Technology and New Needs in recent years, such as Sodium Batteries, Solid-state Battery and Semi-solid-state Battery Equipment, special process requirements such as dry process, 4680 large cylindrical battery equipment requirements etc. Create One-Stop Service for Battery Research Equipment for customers and provide Professional Technical Guidance and Support.

MRBEST related products include: Vacuum Mixers, Planetary Vacuum Mixers, Lab Transfer Coating Machine, Automatic Film Coating Machine, Electric Vertical Roller Machine, Pneumatic Die-cutting Machines, Electric Sealing Machines, Vacuum Pre-sealing machine, Secondary Vacuum Sealing Machine, Coin Cell battery Sealing Machine, Battery Winding Machine, Pouch Cell Sealing Machine, Ultrasonic Spot Welding Machine, Glove Box with Purification System etc. MRBEST products have been sold at home and abroad. Our main customer include: Universities, Research Institutes, Battery Positive and Negative Material Companies, Battery Cell Manufacturers, and New Energy EV Battery Manufacturers etc.





**Shenzhen Meirui Zhida Technology Co., Ltd.**



+86 19866156608



tina@szmrbest.com



lithiumbattery-machine.com

104 Building 6, Second Industrial Zone, Shanmen Community, Yanluo Street, Baoan District, Shenzhen China