




## CR2025 CR2016 CR2032 Coin Cell Cases With Spacer And Spring

Our Product Introduction

for more products please visit us on [lithiumbatterymachine.com](http://lithiumbatterymachine.com)

### Basic Information

- Place of Origin: China
- Brand Name: MRBEST
- Model Number: CR2032/ CR2025/ CR2016
- Minimum Order Quantity: 100 Sets
- Price: Negotiable
- Packaging Details: Standard Export Carton Packing
- Delivery Time: 3-15 working days
- Payment Terms: T/T, Western Union, MoneyGram
- Supply Ability: 5000000 sets per month

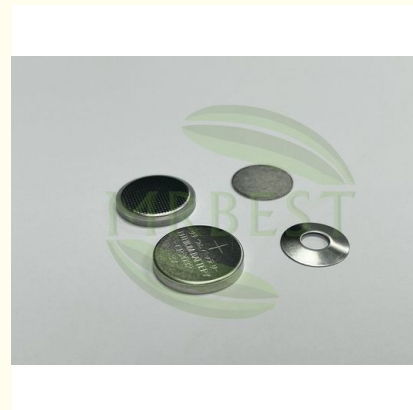


### Product Specification

- Name: Coin Cell Cases
- Shape: Button
- Material: SS304/ SS316
- Battery Size: 20 Mm \*3.2 Mm
- Sealing Gasket: PP
- Full Set Include: Positive Case, Negative Case, Spring And Spacer
- Highlight: **CR2032 Coin Cell Cases, CR2025 Coin Cell Cases**



### More Images

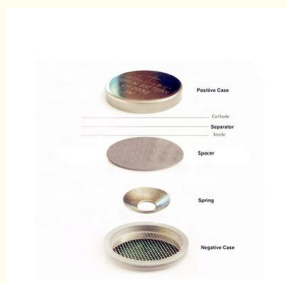

## Product Description






### CR2032 CR2025 CR2016 Coin Cell Cases with Spacer and Spring For Coin Cell Research

CR2032/ CR2025/ CR2016 button cells cases, made of stainless steel 304 or 316 with PP sealing O-ring

Dimensions (Dia x Height)	20 mm x 3.2 mm
Weight	1KG
Material	Casing: SS304or SS316, please choose according to your need. Sealing Gasket: PP The thickness of the SS304/316 for both anode and cathode case is 0.25 mm +/- 0.02 mm.
Application Note	After perfectly crimping, the External Height should be around 3.0mm. The diameter of the inner ring of the plastic gasket for the coin cell 2032 is around 16.39mm. Usually, a case does not need to clean before crimping, but if you do need to clean the case, please put in ethanol for 4-6 hrs, and use Oven to dry several hours, (please set the oven to 30C, if it has ORING.)
Model Choose	CR2032/ CR2025/ CR2016

### Picture Display

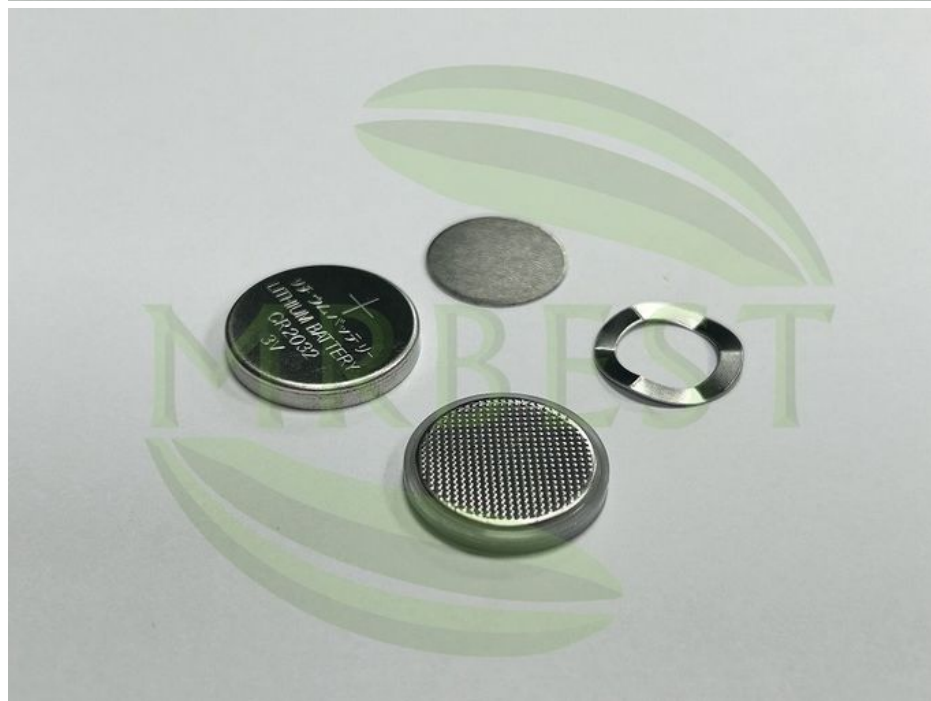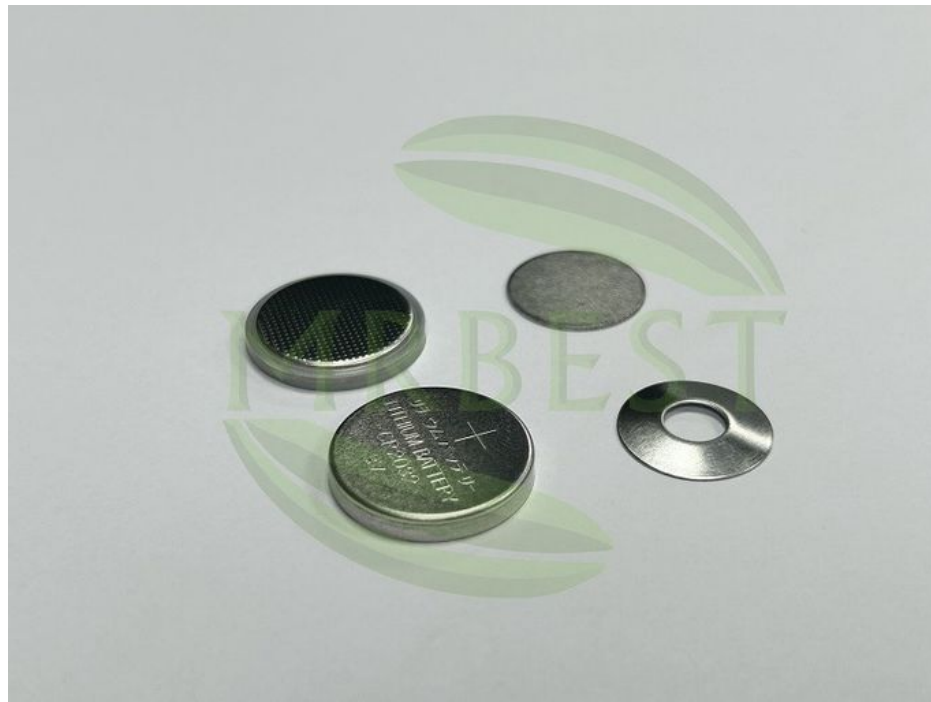

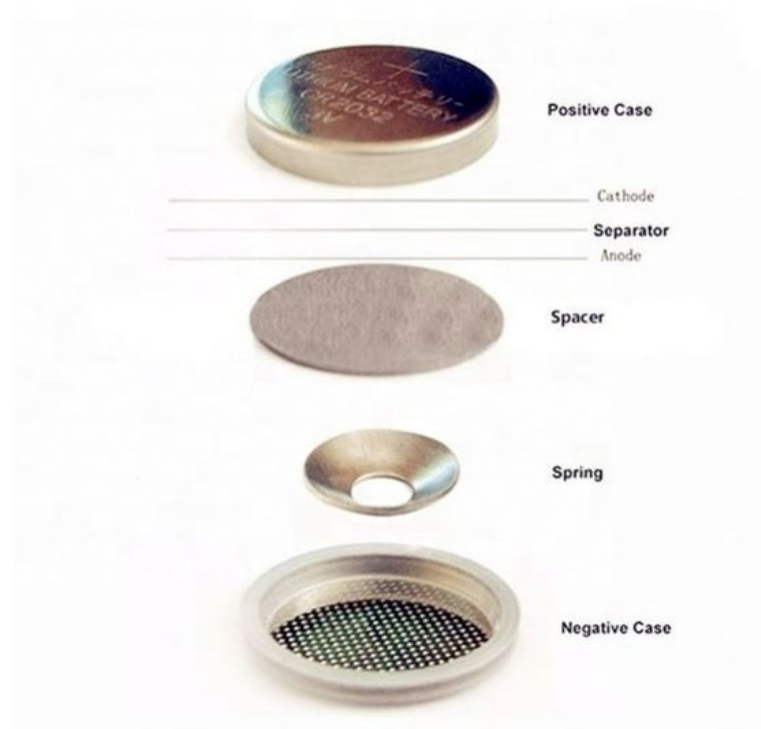











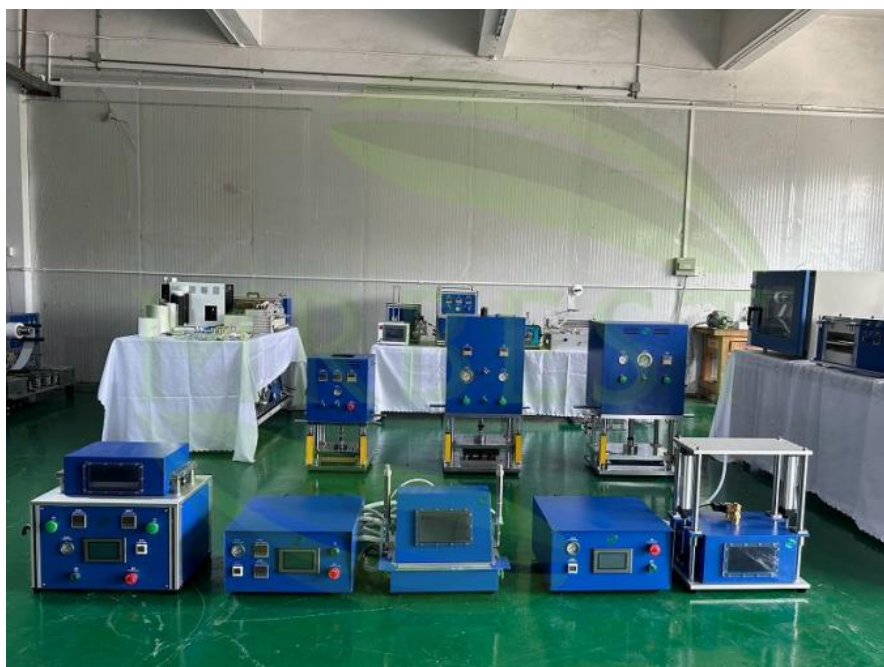

Planetary Vacuum Mixer

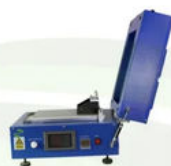

Automatic Film Coating Machine



Rolling Press Machine



Coin Cell Punching Machine

### Coin Cell Lab Equipment



Vacuum Drying Oven



Glove Box



Coin Cell Crimper



Battery Testing Machine





Factory Show



# COMPANY PROFILE

## 20 Years Lithium Battery Industry Experience



Shenzhen Meirui Zhida Technology Co., Ltd. is a Manufacturer specializing in a Complete Set of Lithium Ion Battery Research Equipment. MRBEST has been committed to providing customers with comprehensive Battery Research and Development Equipment to meet customers' new Lithium Battery Technology and New Needs in recent years, such as Sodium Batteries, Solid-state Battery and Semi-solid-state Battery Equipment, special process requirements such as dry process, 4680 large cylindrical battery equipment requirements etc. Create One-Stop Service for Battery Research Equipment for customers and provide Professional Technical Guidance and Support.

MRBEST related products include: Vacuum Mixers, Planetary Vacuum Mixers, Lab Transfer Coating Machine, Automatic Film Coating Machine, Electric Vertical Roller Machine, Pneumatic Die-cutting Machines, Electric Sealing Machines, Vacuum Pre-sealing machine, Secondary Vacuum Sealing Machine, Coin Cell battery Sealing Machine, Battery Winding Machine, Pouch Cell Sealing Machine, Ultrasonic Spot Welding Machine, Glove Box with Purification System etc. MRBEST products have been sold at home and abroad. Our main customer include: Universities, Research Institutes, Battery Positive and Negative Material Companies, Battery Cell Manufacturers, and New Energy EV Battery Manufacturers etc.



**Shenzhen Meirui Zhida Technology Co., Ltd.**



+86 19866156608



tina@szmrbest.com



lithiumbattery-machine.com

104 Building 6, Second Industrial Zone, Shanmen Community, Yanluo Street, Baoan District, Shenzhen China