



Customized Coin Cell Manufacturing Equipment Lithium Ion Battery Making Machine

Our Product Introduction

Basic Information

- Place of Origin: China
- Brand Name: MRBEST
- Certification: CE
- Model Number: MRBEST-Coin Cell
- Minimum Order Quantity: 1 Set
- Price: Negotiable
- Packaging Details: Standard Export Wooden Packing
- Delivery Time: 7-20 working days
- Payment Terms: T/T, Western Union, MoneyGram
- Supply Ability: 500 sets per month



Product Specification

- Product Name: Coin Cell Manufacturing Equipment For Lithium-ion Batteries
- Application: For Battery Research
- Suitable Battery: Coin Cell
- Size: Accepted Customized Size
- Manufacturer: YES
- Direct Factory: With Competitive Price
- High Light: Battery Manufacturing Machine, Battery Making Machine Price, Lithium Battery Making Machine
- Highlight: **Customized Coin Cell Manufacturing Equipment , Coin Cell Lithium Ion Battery Making Machine, Customized Lithium Ion Battery Making Machine**

for more products please visit us on lithiumbatterymachine.com

Product Description

Coin Cell Manufacturing Equipment for Lithium-ion Batteries


Coin/Button Battery Preparation And Assembly

1. Positive Material Ingredients: Active Substance : PVDF : Conductive Agent (Acetylene Black) = 8:1:1
2. Negative Material Ingredients: Active Substance : PVDF : Conductive Agent (Acetylene Black) = 7:2:1
3. Solvent NMP is determined according to the material conditions
4. Coating, use film applicator to do coating
5. Dry, Dry naturally 60 , Vacuum Drying 110 (24h), If water molecules are within the crystal lattice of the material, the drying temperature is 60-70
6. Put the positive and negative electrode pieces into the vacuum glove box for assembly

Notes:

Regardless Of The Positive Electrode Material Or The Negative Electrode Material, They Are Assembled As The Positive Electrode, And The Lithium Chip Is The Negative Electrode.

The Schematic Diagram Of 2016 And 2025 Model Lithium Coin Cell Battery Assembly Is As Follows:

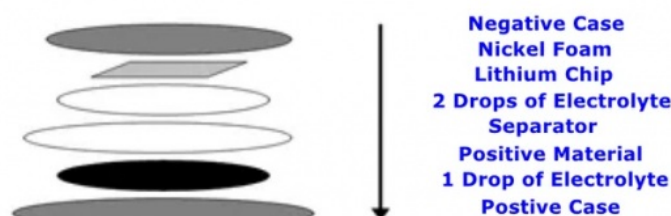

Coin Cell Preparation



MRBEST
 Shenzhen Meirui Zhida Technology Co., Ltd
 Tel/Fax: +86 755 23696396 meiruitop@gmail.com www.mrbest.en.alibaba.com

Our Advantages

1. Steady and high quality ! -cause we are a manufacturer of Lithium Battery Equipment, and we only focus on one thing that is make professional Lithium ion Battery Equipment for customer around the world.
2. Friendly use! all MRBEST Battery Equipment with teach manual and support live teaching via media.
3. Reliable factory! My factory located in Shenzhen city of China, support factory visit and video factory visit .
4. 1 year warranty plus lifetime technical support.
5. All the parts use on machine, in stock forever.
6. Powerful technical support, NOT TRADE COMPANY! We are professional on Lithium Battery Equipment, the chief engineer team have

more than 20 years experience on lithium battery area, so we can easy do diagnose on machine if customer have question.

7. 100% PATIENCE to all customers! Our patience not only before consult order, also after-sale service, one by one .



Shenzhen Meirui Zhida Technology Co., Ltd.



+86 19866156608



tina@szmrbest.com



lithiumbatterymachine.com

104 Building 6, Second Industrial Zone, Shanmen Community, Yanluo Street, Baoan District, Shenzhen China