



Laboratory Doctor Blade Film Coater Coating Machine 250mm Maximum Length

Our Product Introduction

Basic Information

- Place of Origin: China
- Brand Name: MRBEST
- Certification: CE
- Model Number: MR-TMH250
- Minimum Order Quantity: 1 Set
- Price: Negotiable
- Packaging Details: Standard Export Wooden Packing
- Delivery Time: 7-20 working days
- Payment Terms: T/T, Western Union, MoneyGram
- Supply Ability: 500 sets per month




Product Specification

- Name: Laboratory Doctor Blade Film Coater Coating Machine
- Doctor Blade: 100mm Fixed Width As Standard, Other Size Optional
- Coating Length: 250mm Maximum
- Max Heating Temperature: 150
- Heating Type: Top Heating
- Vacuum Pump: Built In Oil Free Vacuum Pump
- Highlight: **Laboratory Doctor Blade Film Coater,
Doctor Blade Film Coater 250mm,
Laboratory Doctor Blade Film Coating Machine**

for more products please visit us on lithiumbatterymachine.com

Product Description

Laboratory Doctor Blade Film Coater Coating Machine

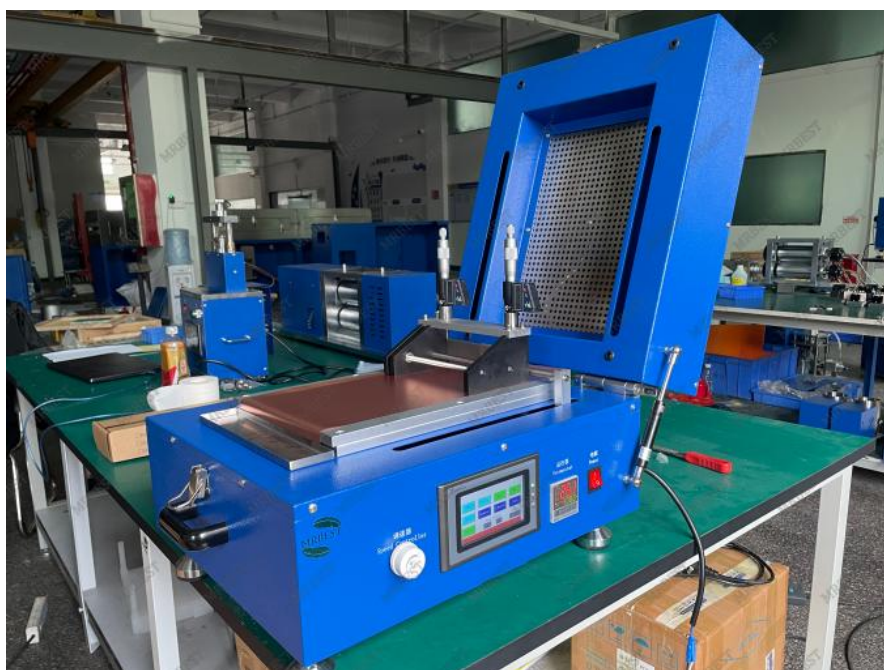

MR-TMH250 is a compact film tape casting machine with vacuum chuck, which widely used in all kinds film coating research, like ceramic films, crystal films, lithium battery electrode coating, etc. It is designed with top cover heating, the heating temperature up to 150 adjustable and digital display.

Features

- Touch screen operation, high automation.
- Vacuum aluminium chuck is designed for quick suction of substrate.
- Built-in vacuum pump is designed.
- Coating stroke length can be set.
- Heating temperature up to 150 adjustable, and it's digital display.

Technical Parameters

Working voltage	AC220V
Coating length	Max 250mm adjustable
Doctor blade	100mm fixed width as standard, other size is optional
Coating thickness	0.01~3.5mm adjustable
Vacuum chuck	L365 x W200 x H30mm
Vacuum pump	Built In Oil Free Vacuum Pump
Coating speed	0~200mm/s adjustable
Heating method	Top cover heating
Heating temperature	RT~150 adjustable, digital display
Heating temperature accuracy	±1
Dimension	L500 x W310 x H330mm
Weight	45kg

Our Advantages

1. Steady and high quality ! -cause we are a manufacturer of Lithium Battery Equipment, and we only focus on one thing that is make professional Lithium ion Battery Equipment for customer around the world.
2. Friendly use! all MRBEST Battery Equipment with teach manual and support live teaching via media.
3. Reliable factory! My factory located in Shenzhen city of China, support factory visit and video factory visit .
4. 1 year warranty plus lifetime technical support.
- 5.All the parts use on machine, in stock forever.
6. Powerful technical support, NOT TRADE COMPANY! We are professional on Lithium Battery Equipment, the chief engineer team have more than 20 years experience on lithium battery area, so we can easy do diagnose on machine if customer have question.
7. 100% PATIENCE to all customers! Our patience not only before consult order, also after-sale service, one by one .





+86 19866156608



tina@szmrbest.com



lithiumbattery-machine.com

104 Building 6, Second Industrial Zone, Shanmen Community, Yanluo Street, Baoan District, Shenzhen China